



## CryoCarrier Series

### Dry Shipper for Transportation

CryoCarrier Series is dry shipper dewars for safe transportation of biological products in vapor phase at cryogenic (-150°C or lower) temperatures.

Taking advantage of innovative absorbing material and optimized construction, Antech CryoCarrier Series dry shippers can be recharged within 2 hours, while new/warm tanks still need a sit time of 24 hours.

### Features

- Vapor phase storage ensures dry and spill-free shipping
- Light in weight and durable to use
- Applicable for various of sample containers, such as vial, straw, blood bag
- Recharged in 2 hours
- Holdover time 2~15 days to meet different requirements
- Protective trolley case as optional

Based on type of products to be transported, Antech offers 3 lines of dry-shippers:

- CryoCarrier R: for biological products stored with cryogenic vials and boxes;
- CryoCarrier B: for biological products stored with 25ml/50ml/250ml blood bags and cassettes;
- CryoCarrier C: for biological products stored with round canister, such as straws, cryogenic vials.



CryoCarrier R



CryoCarrier B



CryoCarrier C

## Protective Trolley Case

- Heavy duty case with aluminum frame, durable to use
- Forming material protection for safety
- Trolley for easy movement
- Case lock preventing unauthorized access




## Smart Data Logger

- -200~70°C
- Multi-Use, replaceable battery
- USB 2.0 plug, 28800 readings
- Encrypted PDF report embedded with raw data Bluetooth printer connection (200XB)



## Specification

CryoCarrier R: for biological products stored with cryogenic vials and boxes;

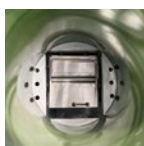


Model	CryoCarrier 6R	CryoCarrier 10R	CryoCarrier 15R	CryoCarrier 25R
<b>Maximum Storage Capacity</b>				
2ml vial	200	100	300	500
25-place box	—	4	—	—
81/100-place box, 2 inch height	2	—	3	5
5ml vial	—	50	81	162
81-place box, 3 inch height	—	—	1	2
<b>Performance</b>				
Container gross capacity, liter	6	10	15	25
LN <sub>2</sub> absorption, liter	1.8	6	6	10
Static evaporation rate (L/day)	0.8	0.4	0.89	0.89
Vapor working duration (day)	2	15	7	11
<b>Dimensions</b>				
Neck Opening (mm)	190	125	216	216
External diameter (mm)	252	300	396	396
Overall height (mm)	519	584	610	741
Weight empty (kg)	7.2	8.5	16.0	22.0
Weight charged, vapor (kg)	8.7	13.4	20.9	31.8
<b>Accessories</b>				
Protection trolley case, optional	DSC-6	DSC-10	DSC-15	DSC-25
Data logger, optional	200X, 200XB			

CryoCarrier B: for biological products stored with 25ml/50ml/250ml/500ml blood bags and cassettes.



Model	CryoCarrier 4B	CryoCarrier 6B	CryoCarrier 25B
<b>Maximum Storage Capacity</b>			
25ml blood bag	6	5	17
50ml blood bag	6	5	11
250ml blood bag	—	5	11
500ml blood bag	—	5	11
<b>Performance</b>			
Container gross capacity, liter	4	6	25
LN <sub>2</sub> absorption, liter	2	1.8	10
Static evaporation rate (L/day)	0.45	0.8	0.89
Vapor working duration (day)	4	5	11
<b>Dimensions</b>			
Neck Opening (mm)	125	190	216
External diameter (mm)	300	252	396
Overall height (mm)	521	519	741
Weight empty (kg)	7.0	7.2	22.0
Weight charged, vapor (kg)	8.6	8.7	31.8
<b>Accessories</b>			
Protection trolley case, optional	DSC-4	DSC-6	DSC-25
Data logger, optional	200X, 200XB		

CryoCarrier C: for biological products stored with round canister, such as straws, cryogenic vials;

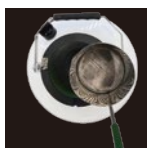

Model	CryoCarrier 2C	CryoCarrier 3C	CryoCarrier 4C	CryoCarrier 6C	CryoCarrier 10C	CryoCarrier 25C
<b>Maximum Storage Capacity</b>						
1/4 cc straw	115	451	1422	1065	2046	11584
1/2 cc straw	53	181	664	492	942	5456
2ml vial	3	21	84	66	132	678
<b>Performance</b>						
Container gross capacity, liter	1.5	3	4	6	10	25
LN <sub>2</sub> absorption, liter	0.6	1.3	1.8	2.4	5.8	7
Static evaporation rate (L/day)	0.07	0.15	0.19	0.21	0.35	0.89
Vapor working duration (day)	9	9	9	11	17	8
<b>Dimensions</b>						
Neck Opening (mm)	30	50	70	80	80	216
External diameter (mm)	180	224	224	300	300	396
Overall height (mm)	398	441	505	492	584	741
Canisters diameter (mm)	25	47	67	76	76	190
Canister height (mm)	120	200	276	200	280	276
Weight empty (kg)	3.0	4.0	5.5	7.0	7.7	22.0
Weight charged, vapor (kg)	3.5	5.1	7.0	9.0	12.4	27.7
<b>Accessories</b>						
Protection trolley case, optional	—	DSC-4	DSC-4	DSC-6	DSC-10	DSC-25
Data logger, optional	200X, 200XB					